



Using Steps

· Step 1: Fix the toothbrush head on the handle correctly.

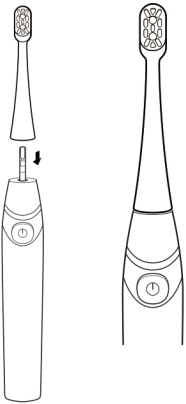


· Step 2: Squeeze proper amount of toothpaste on toothbrush head.

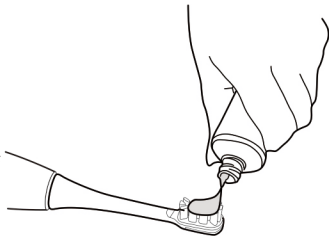

5

Steps 3

· Step 3: Move slowly to brush different areas of your oral cavity; press the multi-functional button to switch between 3 modes; short frequency vibration during brushing reminds you to change cleaning area.


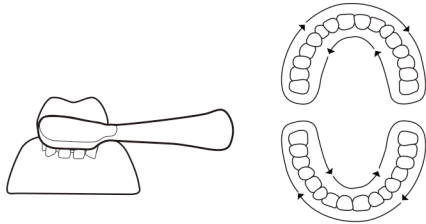


Note:

- To avoid toothpaste falling off, put toothbrush head into your mouth before turning on the toothbrush.
- To avoid foam splashing, turn off the toothbrush before taking toothbrush head out of your mouth.

6


After cleaning your teeth, wash the toothbrush head with flowing water. Cover the toothbrush head with the toothbrush cap.

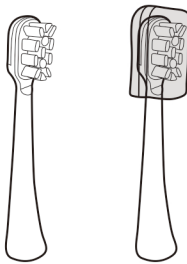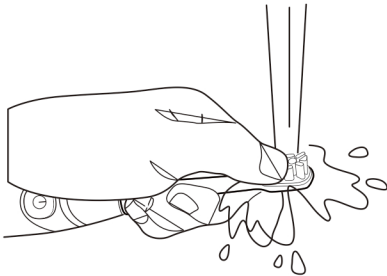

7

Toothbrush Sterilization

To keep your toothbrush hygienic, we advise to sterilize your toothbrush head at least once a week.
· Take off the toothbrush cap and keep the water out before sterilization. Insert the toothbrush head into the hole. The sterilization and drying function will be automatically on after 2 seconds.
· The purple breathing light will be on when the sterilization and drying function work. Sterilization function will turn off automatically after 15 minutes. Then the indicator light will turn blue. Drying function needs 15 minutes.

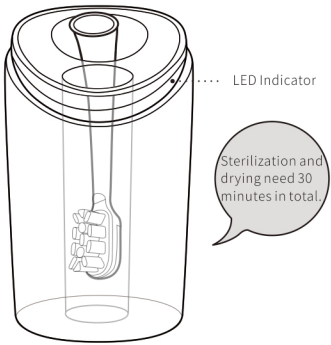


8

Product Specifications

Main body
Main body name: Allegro Sonic Electric Toothbrush
Model: al-01
In-put: 5V=1A
Power: 1W±0.7W
Water-proof level: IPX7
Origin: China

Multi-functional charging kit
Main function: Inductive charging/sterilization/drying
Rated in-put: 5V-1A
Sterilization time: 15 minutes
Drying time: 30minutes(including sterilization time)
Water-proof level: Not water-proof
Charging time: 16 hours to be fully charged

FC RoHS CE

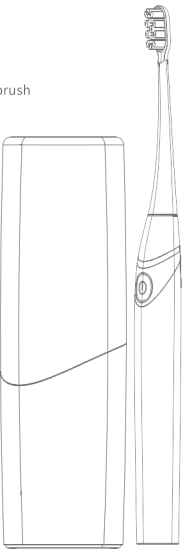

9

Aellgro

Sonic Electric Toothbrush

User Manual

Caution:

1. There could be gingival tingling, bleeding and other symptoms at first use. Please use sensitive mode for 1 to 2 weeks. If you still feel uncomfortable after 2 weeks, there could be something wrong with your gums. We advise you to stop using the product and consult your dentist about the problem.
2. Use the product under your dentist's guidance if you have had an orthodontic treatment, dental repairing or dental implants etc.
3. Don't use the product for other purpose other than cleaning your teeth.
4. Make sure your toothbrush and charging base won't fall off from where you put them. If they fall down, please check if they are still water-proof.
5. Don't disassemble or repair the product by yourself. It may expose you in danger.

Prohibition


1. Don't put the charging base in water or any other liquid.
2. Please confirm that plug is dry when charging and sterilizing. Don't use the product in humid environment. Turn off the power if you are not using it for a long time.
3. This product contains a built-in rechargeable battery. Don't put, use or charge it in high temperature.
4. If the product gets damaged, please stop using it.
5. Keep the product away from strong magnetic field (such as induction cooker) when charging.
6. Please don't stare at the purple indicator light which may hurt your eyes.

Reminder

1. You can share the toothbrush main body with others but not the toothbrush head.
2. We recommend you to replace the toothbrush head once every three months.
3. To guarantee your safety and experience, please use accessories from Allegro store.
4. To extend the battery's life of the toothbrush, please get your toothbrush fully discharged and then charge it for more than 16 hours once every 6 months.

1


Product description

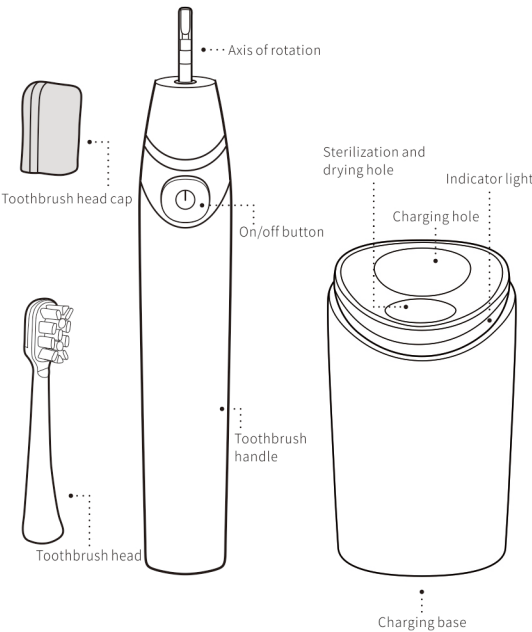

2

First Using and Charging

- Please charge the electric toothbrush for at least 16 hours before first use.
- When charging, a red breathing light is on. When fully charged, the indicator light turns blue.
- When power on, there is a blue indicator light on the charging base.
- The sterilization and drying functions have an inductive switch. No manual operation is needed.
- When sterilization function works, a purple breathing light is on. Sterilization function will automatically turn off after 15 minutes and drying function will automatically turn off after 30 minutes.

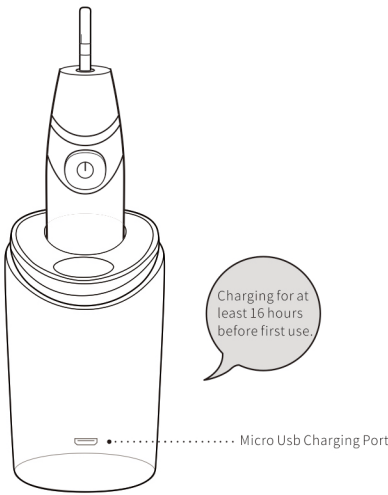


3

Choosing Mode

Turn on the Allegro Electric Toothbrush and switch between 3 modes with the multi-functional button.
Press once: clean mode, the indicator light will flash at a moderate speed.
Press two times: sensitive mode, the indicator light will flash slowly.
Press three times: clean and white mode, the indicator light will flash quickly.
Press four times: turn off the toothbrush.

Note:

- Choosing a mode should be done within 20 seconds after turning on the toothbrush.
- Pressing on the multi-functional button beyond 20 seconds means turning off the toothbrush.
- Press on the multi-functional button for 2 seconds to turn off the toothbrush in any mode.

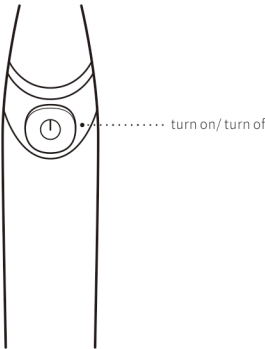

4